


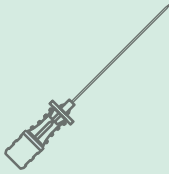
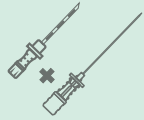
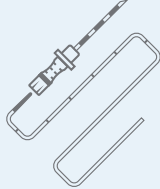
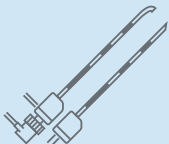




NRFit™ Brochure

OUR VISION

We
PROTECT
and
IMPROVE
the
HEALTH
of people
around the world.

CONTENT

<p>What is NRFit™?</p>		<p>WHAT IS NRFit™?</p>
<p>SPINAL NEEDLES AND PACKS</p> <p>Spinocan® Pencan® Spinal Introducer Needle</p>		<p>SPINAL NEEDLES AND PACKS</p>
<p>COMBINED SPINAL EPIDURAL</p>		<p>COMBINED SPINAL EPIDURAL</p>
<p>EPIDURALS AND PACKS</p> <p>Perican® Perican® Paed Epican® Paed Perifix® Epidural Catheter Accessories</p>		<p>EPIDURALS AND PACKS</p>
<p>PERIPHERAL NERVE BLOCKS</p> <p>Stimuplex® Ultra 360® Ultraplex® 360 Contiplex® C Contiplex® S/S Ultra 360®</p>		<p>PERIPHERAL NERVE BLOCKS</p>
<p>NRFit™ ACCESSORIES</p>		<p>NRFit™ ACCESSORIES</p>
<p>TRAINING AND SUPPORT</p>		<p>TRAINING AND SUPPORT</p>



1. WHAT IS NRFit™?

THE REASON

- The universally compatible Luer connector has been identified in wrong route administration.
- There have been a number of errors recorded where patients have been severely harmed or died because of simple misconnections.¹

THE RESPONSE


The National Patient Safety Agency published a series of alerts in 2009 and 2011 requiring NHS trusts to use non-luer devices for neuraxial procedures.

THE SOLUTION

The International Organisation for Standardisation (ISO) has developed a standard for small bore connectors in Neuraxial and Regional Anaesthesia (ISO 80369-6), called NRFit™.

HOW DOES NRFit™ DIFFER?

- Connectors are 20% smaller in diameter.
- They will only connect to other NRFit™ Devices.
- They are yellow in colour where feasible.

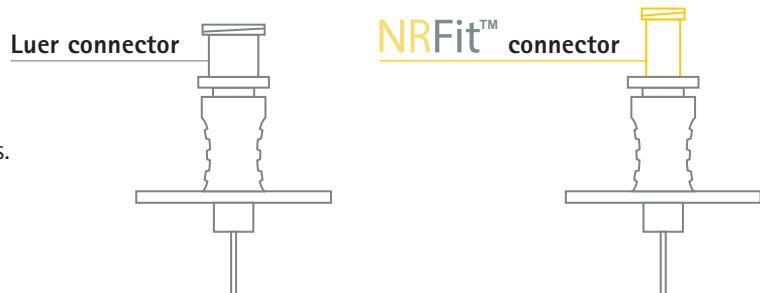

WHAT IS NRFit™?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL


EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFit™ ACCESSORIES

TRAINING AND SUPPORT

SPINAL NEEDLES AND PACKS

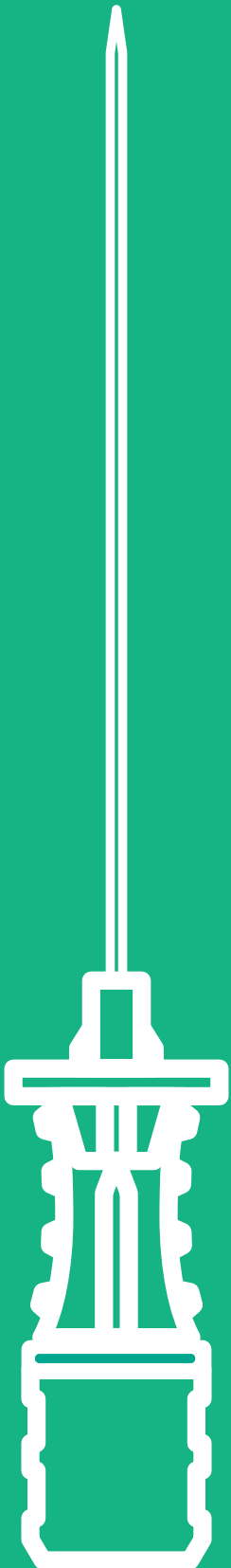

2. SPINAL NEEDLES AND PACKS


PENCAN®

Atraumatic tip design results in low PDPH rates and superior tactile feedback.

Shorter distance from tip to start of eyelet (1.05 mm), designed for quicker flash back.


Shorter distance from tip to back of eyelet (1.7 mm), designed to minimise insertion depth of the needle tip into the subarachnoid space.

Eyelet is sized at 0.65 mm to help reduce risk of inaccurate drug delivery caused by 'straddling the dura'.

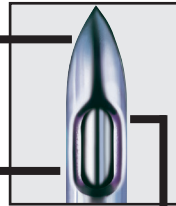

Colour coded hub to identify gauge size.

INTRODUCERS

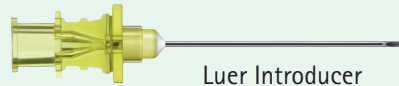

Luer Introducer

B. Braun Code	4505000
NHS Supply Chain Code	FTR2707

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Pencan 22G x 88 mm	4502035N-01	FTR2368	
NRFit Pencan 25G x 25 mm	4502167N-01	FTR2599	
NRFit Pencan 25G x 50 mm	4502159N-01	FTR2375	
NRFit Pencan 25G x 88 mm	4502019N-01	FTR2600	
NRFit Pencan 25G x 88 mm with Luer Introducer	4502043N-01	FTR2780	
NRFit Pencan 25G x 103 mm with Luer Introducer	4502116N-01	FTR2781	
NRFit Pencan 25G x 120 mm	4502120N-01	FTR2369	25
NRFit Pencan 25G x 156 mm	333877N-01	FTR2370	
NRFit Pencan 27G x 50 mm	4502175N-01	FTR2601	
NRFit Pencan 27G x 88 mm	4502027N-01	FTR2699	
NRFit Pencan 27G x 88 mm with Luer Introducer	4502051N-01	FTR2710	
NRFit Pencan 27G x 103 mm with Luer Introducer	4502124N-01	FTR2711	
NRFit Pencan 27G x 120 mm	4502132N-01	FTR2371	

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT

2. SPINAL NEEDLES AND PACKS

SPINOCAN®



Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Spinocan 18G x 75 mm	4501373N-01	FTR2585	
NRFit Spinocan 18G x 88 mm	4501390N-01	FTR2360	
NRFit Spinocan 19G x 40 mm	4501144N-01	FTR2586	
NRFit Spinocan 19G x 88 mm	4501195N-01	FTR2587	
NRFit Spinocan 20G x 75 mm	4509757N-01	FTR2588	
NRFit Spinocan 20G x 88 mm	4509900N-01	FTR2361	
NRFit Spinocan 22G x 40 mm	4507401N-01	FTR2362	
NRFit Spinocan 22G x 75 mm	4507754N-01	FTR2363	
NRFit Spinocan 22G x 88 mm	4507908N-01	FTR2364	
NRFit Spinocan 22G x 120 mm	4506090N-01	FTR2365	
NRFit Spinocan 22G x 150 mm	4506095N-01	FTR2366	
NRFit Spinocan 25G x 75 mm	4505751N-01	FTR2590	
NRFit Spinocan 25G x 88 mm	4505905N-01	FTR2367	
NRFit Spinocan 25G x 120 mm	4505913N-01	FTR2589	
NRFit Spinocan 26G x 88 mm	4502906N-01	FTR2592	
NRFit Spinocan 26G x 120 mm	4504917N-01	FTR2591	
NRFit Spinocan 27G x 88 mm	4503902N-01	FTR2594	
NRFit Spinocan 27G x 120 mm	4502140N-01	FTR2593	
NRFit Spinocan 29G x 88 mm with Luer Introducer	4501900N-01	FTR2772	
NRFit Spinocan 29G x 120 mm	4501918N-01	FTR2595	

NRFit™ BASIC PROCEDURE PACK



LUER CONTENTS

- Sterican Hypodermic Needle 18G x 40 mm
- Sterican Hypodermic Needle 25G x 16 mm
- Omnifix® 5 ml Slip Syringe

NRFit™ CONTENTS

- NRFit™ Filter Needle 5µm
- NRFit™ Omnifix® 1 ml Slip Syringe
- NRFit™ Omnifix® 3 ml Slip Syringe
- NRFit™ Omnifix® 5 ml Slip Syringe

ACCESSORIES

- Transparent Drape (75 cm x 11 cm)
- Gauze Compresses x 5
- Sterile Field Sheet (90 cm x 90 x cm)
- Green Tray

B. Braun Code	8864006A
NHS SC Code	FTR2839
Unit of Measure	12

NRFit™ SPINAL MINI PACK



CONTENTS

- NRFit™ Pencan® 25G x 88 mm Spinal Needle with Luer Introducer
- NRFit™ Filter Needle 5µm
- NRFit™ Omnifix® 5 ml Slip Syringe

B. Braun Code	4470000N
NHS SC Code	FTR3292
Unit of Measure	100

WHAT IS
NRFIT™?

SPINAL NEEDLES AND
PACKS

COMBINED
SPINAL EPIDURAL


EPIDURALS AND PACKS

PERIPHERAL NERVE
BLOCKS

NRFIT™
ACCESSORIES

TRAINING AND
SUPPORT

COMBINED SPINAL EPIDURAL

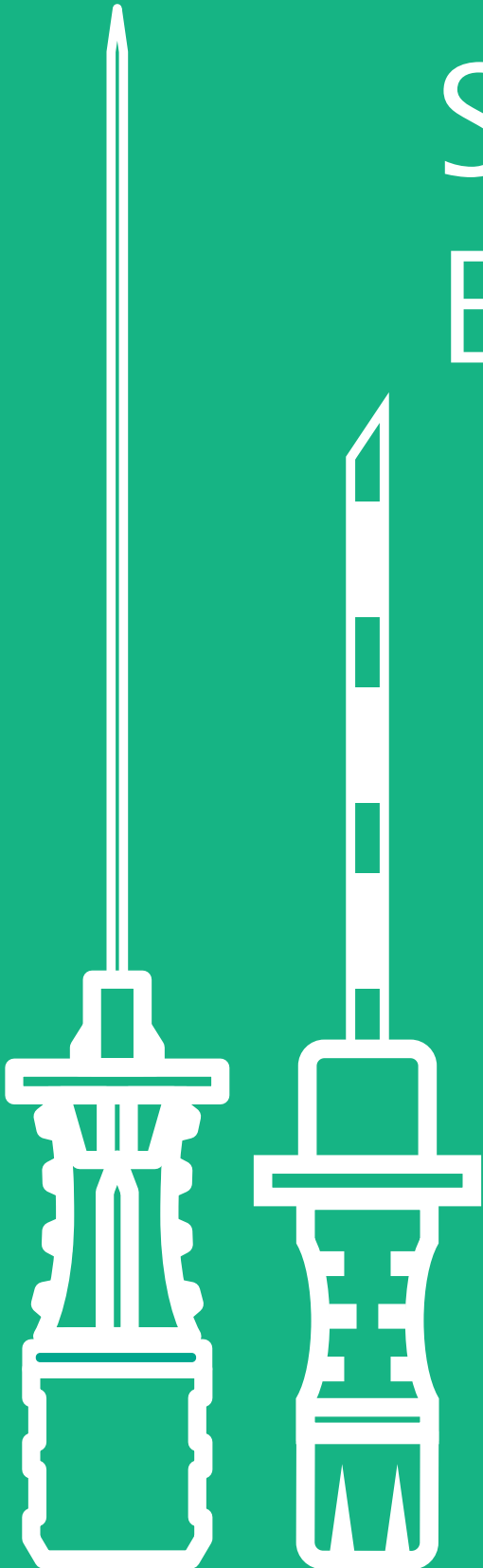

3. COMBINED SPINAL EPIDURAL



NRFit™ CONTENTS

- NRFit™ Perican® with Back Eye
- NRFit™ Pencan® or Spinocan®
- NRFit™ Catheter + Thread Assist Guide
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix Pin Pad
- NRFit™ Omnifix® Loss of Resistance 8 ml Syringe

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Soft-Tip with Docking Station + Pencan 27G x 141 mm	4556763N-01	FYB1098	10
NRFit Soft-Tip with Docking Station + Spinocan 27G x 139 mm	4556747N-01	FYB1095	
NRFit Standard Espocan with Pencan 27G x 130.5 mm	4556666N-01	FYB1097	
NRFit Standard Espocan with Spinocan 27G x 126 mm	4556674N-01	FYB1096	

The ideal overlap of the spinal needle is by

Pencil Point Needle


|— 15 mm —|

"Quincke" Needle

|— 13 mm —|

DOCKING SYSTEM

- Ensures fixation of the spinal needle inside the Tuohy needle after successful dural puncture.
- Allows rotation of the spinal needle even in the fixed position.
- Comfortable finger grip to steady spinal needle in position.

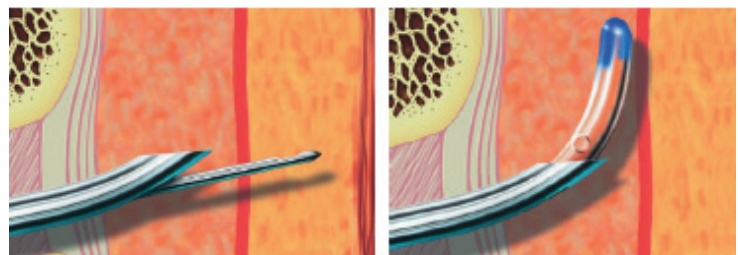

WHAT IS NRFit™?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL


EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFit™ ACCESSORIES


TRAINING AND SUPPORT

EPIDURALS AND PACKS

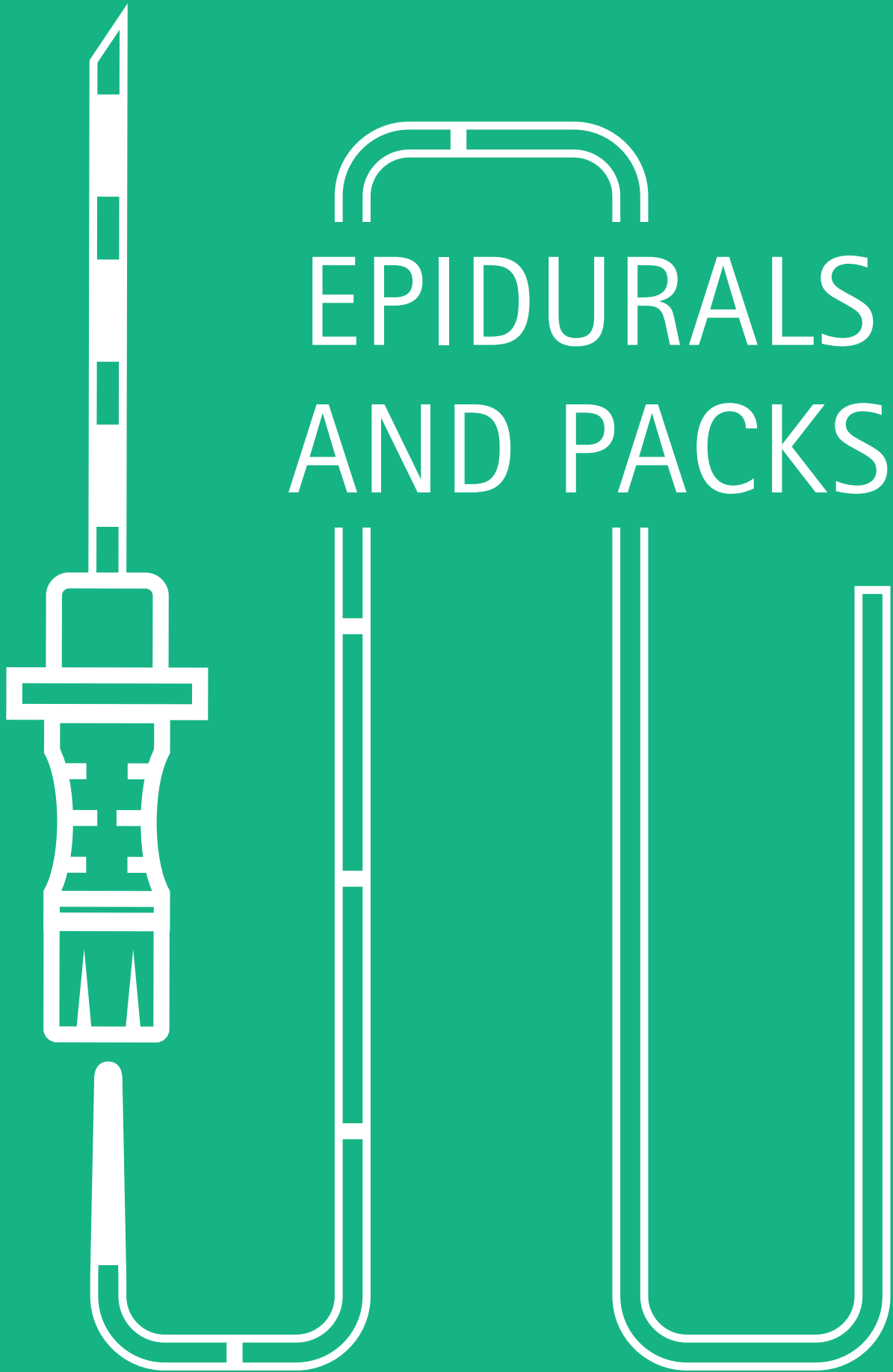

4. EPIDURALS AND PACKS

PERICAN[®] NEEDLE

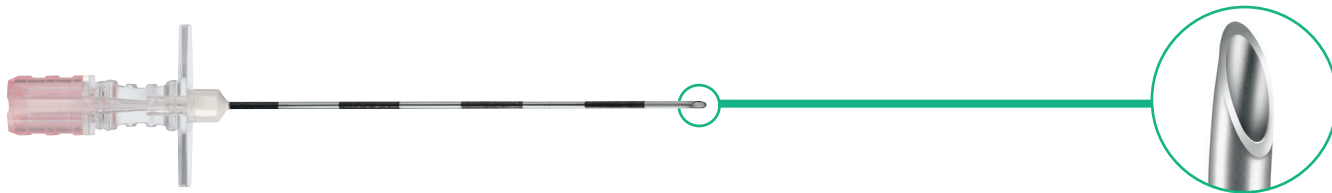

- Polished inner bevel helps to minimise the risk of catheter shearing.
- 1 cm graduation markings all the way down the needle shaft.
- Fixed wings that give extra stability.
- Threading assist guide fits securely into the needle hub and assists easy advancement of the catheter.

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perican 16G x 80 mm	4512782N-01	FTR2462	
NRFit Perican 16G x 120 mm	4512785N-01	FTR2463	
NRFit Perican 16G x 150 mm	4512784N-01	FTR2464	
NRFit Perican 18G x 80 mm	4512383N-01	FTR2465	
NRFit Perican 18G x 120 mm	4512453N-01	FTR2466	25
NRFit Perican 18G x 150 mm	4512200N-01	FTR2467	
NRFit Perican Paed 18G x 50 mm	4502302N-01	FTR2582	
NRFit Perican Paed 20G x 50 mm	4502094N-01	FTR2583	
NRFit Perican Paed 22G x 50 mm	4502078N-01	FTR2584	

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT[™] ACCESSORIES

TRAINING AND SUPPORT

4. EPIDURALS AND PACKS

CATHETER TYPES



THE SOFT TIP CATHETER

- A moderately stiff but flexible shaft with a soft atraumatic tip
- 20G x 1000 mm
- 3 lateral eyes
- Thread assist guide

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix Soft Tip 20G (side openings)	4515048N-01	FTR2774	25



THE STANDARD CATHETER

- Moderately stiff catheter shaft for easy catheter insertion
- 20G or 19G x 1000 mm
- 3 lateral eyes or 1 central opening
- Thread assist guide

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix Standard Catheter 19G (side opening)	4513258N-01	FTR2777	25
NRFit Perifix Standard Catheter 20G (side opening)	4513150N-01	FTR2778	25
NRFit Perifix Standard Catheter 20G (central opening)	4513177N-01	FTR2779	25



THE ONE CATHETER

- Stiff catheter for easy catheter insertion
- 6 lateral eyes help to enhance the spread of the anaesthetic drug
- 3 pairs of micro-holes
- 19G or 20G x 1000 mm
- Thread assist guide

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix One 19G (side openings)	45132581N-01	FTR2775	25
NRFit Perifix One 20G (side openings)	45131501N-01	FTR2776	25

NRFit™ PERIFIX FILTER AND COMPLETE SETS



NRFit™ CONTENTS

- NRFit™ Perican®
- NRFit™ Catheter + Thread Assist Guide
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- NRFit™ Omnifix® Loss of Resistance 8 ml Slip Syringe
- Epidural Sticker

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix One 401 18G Filter Set	45140171N-01	FTR2409	10
NRFit Perifix One 402 16G Filter Set	45140251N-01	FTR2410	
NRFit Perifix Soft Tip 701 18G Filter Set	4510097N-01	FTR2413	
NRFit Perifix Soft Tip 730 16G Filter Set	4517309N-01	FTR2414	
NRFit Perifix Standard 16G 402 Filter Set	4514025N-01	FTR2412	
NRFit Perifix Standard 18G 400 Filter Set (No LOR Syringe)	4514009N-01	FTR2700	
NRFit Perifix Standard 18G 401 Filter Set	4514017N-01	FTR2411	



NRFit™ CONTENTS

- NRFit™ Perican® 18G x 80 mm
- NRFit™ Catheter + Thread Assist Guide
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- NRFit™ Filter Needle 5µm
- NRFit™ Omnifix® 20 ml Lock Syringe
- NRFit™ Omnifix® Loss of Resistance 8 ml Slip Syringe
- Epidural Sticker

LUER CONTENTS

- Omnifix® 3 ml Lock Slip Syringe
- Sterican® 25G x 16 mm Hypodermic Needle
- Sterican® 18G x 40 mm Hypodermic Needle
- Sterican® 21G x 40 mm Hypodermic Needle

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix One 421 18G Complete Set	45142111N-01	FTR2704	10
NRFit Perifix Soft Tip 901 18G Complete Set	4510305N-01	FTR2702	
NRFit Perifix Standard 18G 421 Complete Set	4514211N-01	FTR2701	

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT

4. EPIDURALS AND PACKS

NRFit™ PERIFIX PROCEDURE PACKS

NRFit™ PERIFIX® 18G SOFT-TIP PROCEDURE PACK



NRFit™ CONTENTS

- NRFit™ Perican® 18G x 80 mm
- NRFit™ Soft-Tip Catheter + Thread Assist Guide
- NRFit™ Filter Needle 5µm
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- NRFit™ Omnifix® 1 ml Slip Syringe
- NRFit™ Omnifix® 10 ml Slip Syringe
- NRFit™ Omnifix® Loss of Resistance 8 ml syringe
- Epidural Sticker

LUER CONTENTS

- Omnifix® 5 ml Luer Slip Syringe
- Sterican® 18G x 40 mm Hypodermic Needle
- Sterican® 25G x 16 mm Hypodermic Needle

ACCESSORIES

- Transparent Drape with 10 cm hole (75 cm x 110 cm)
- Hand Towel
- Gauze x 5
- Swab Sticks x 3
- Sterile Field Sheet (90 cm x 90 cm)
- Green Tray

Product Description

NRFit Perifix Soft Tip 18G Procedure Pack

B. Braun Code

8864003

NHS SC Code

FTR3291

UOM

8

NRFit™ PERIFIX® 16G SOFT-TIP PROCEDURE PACK



NRFit™ CONTENTS

- NRFit™ Perican® 16G x 80 mm
- NRFit™ Soft-Tip Catheter 20G
- NRFit™ Filter Needle 5µm
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2 µm
- Perifix® Pin Pad
- NRFit™ Perifix Loss of Resistance Syringe 8 ml Slip
- NRFit™ Omnifix® 1 ml Slip
- NRFit™ Omnifix® 5 ml Lock
- NRFit™ Omnifix® 10 ml Slip

LUER CONTENTS

- Luer Omnifix® 5 ml Slip Syringe
- Sterican® 18G x 40 mm Hypodermic Needle
- Sterican® 25G x 16 mm Hypodermic Needle

ACCESSORIES

- Transparent Drape with 10 cm hole (75 cm x 110 cm)
- Gauze x 5
- Sterile Field Sheets (100 cm x 100 cm)
- Fenestrated Drape (70 cm x 110 cm)
- Green Tray

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perican16G Procdedure Pack	8864003N	TBC	8

NRFit™ PAEDIATRIC SETS



NRFit™ CONTENTS

- NRFit™ Paediatric Perican®
- NRFit™ ONE Catheter + Thread Assist Guide
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- NRFit™ Omnifix® Loss of Resistance syringe
- NRFit™ Omnifix® 3 ml Lock Syringe
- NRFit™ Filter Needle 5µm
- Epidural Sticker

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix One Paed 18G Set	45120061N-01	FTR2705	10
NRFit Perifix One Paed 20G Set	45120141N-01	FTR2706	10

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT

4. EPIDURALS AND PACKS


EPIDURAL ACCESSORIES



LOSS-OF-RESISTANCE (LOR) SYRINGES

- Low friction plunger for smooth glide characteristics
- Parabolic graduation markings to distinguish from conventional syringes

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit LOR 8 ml Slip Syringe	4637110	FWC537	100
NRFit LOR 8 ml Lock Syringe	4638110	FWC536	100

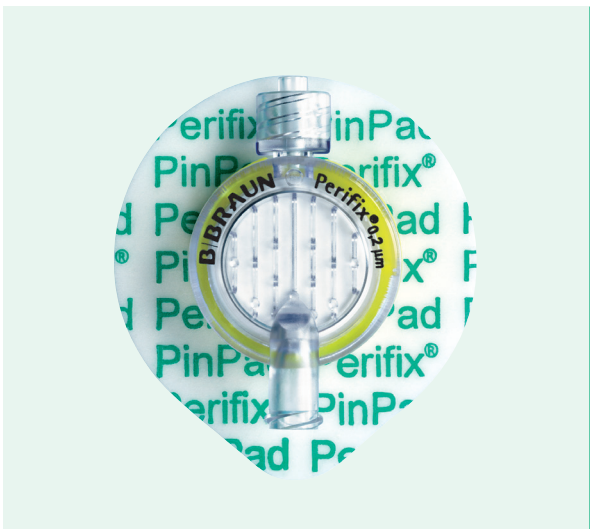

PERIFIX FILTER WITH PIN PAD

- 0.2µm particulate filter
- Flat design for patient comfort
- Self-adhesive pin pad
- Rotatable feature allows the patient to move freely with reduced risk of catheter disconnection


Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix Filter 0.2µm	4515501N-01	FTR2703	25
Perifix Pin pad Unsterile	4515510	FYB861	25



CATHETER CONNECTOR

- Securely and easily fixes catheter in place with a simple snap action
- Audible click verifies solid connection has been made

Product Description	B. Braun Code	UOM
Perifix Catheter Connector NRFit 19G	4513801N-01	25
Perifix Catheter Connector NRFit 20/24G	4513800N-01	25

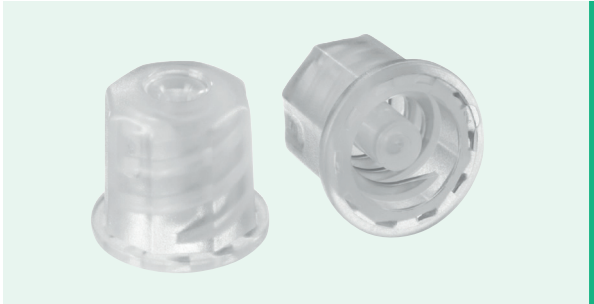

PERIFIX STOPPER

- NRFit™ stoppers are there to close our NRFit™ filter or perifix catheter connector

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Perifix Stopper	4517502N	FWC563	25

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL


EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT

PERIPHERAL NERVE BLOCKS

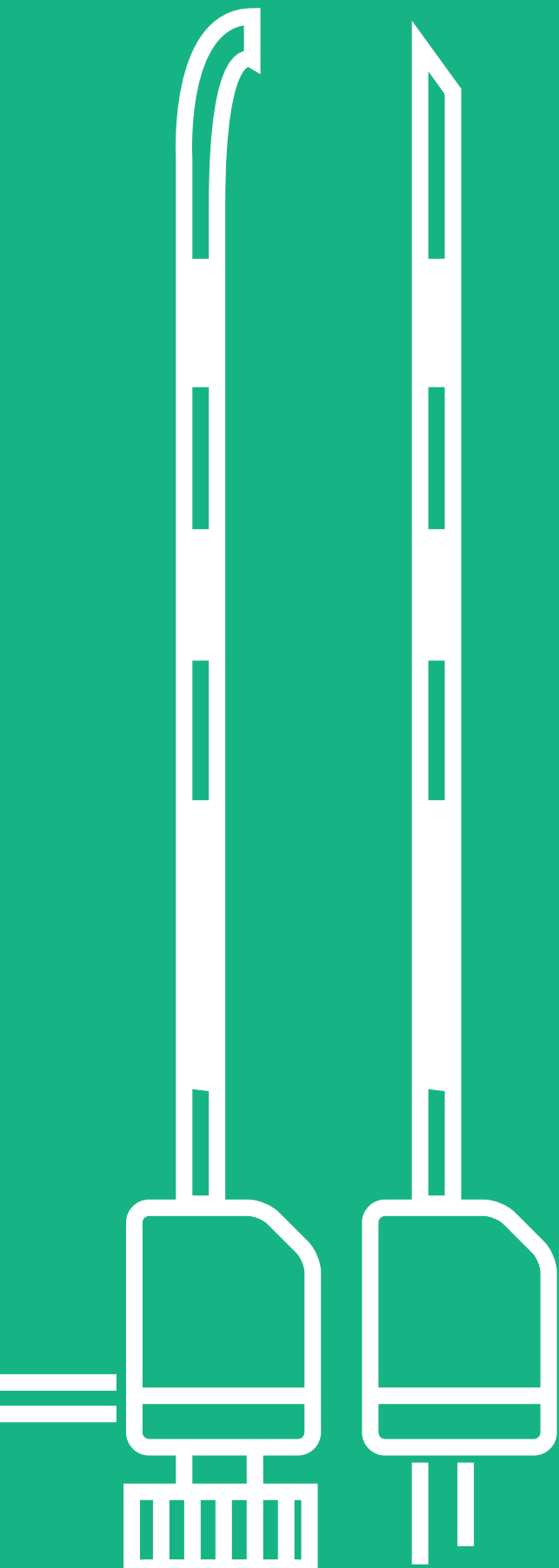

5. PERIPHERAL NERVE BLOCKS

SINGLE SHOT

30 degree bevel for optimal puncture and gliding behaviour.

Laser-crafted 360 echogenic X-pattern facilitates clear needle visualisation from all sides.

Safety code enables easier identification of the needle and needle tip position.

Stimuplex® (Stimulating) or Ultraplex® (Non-Stimulating).

50 cm long DEHP-free injection tubing.

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Stimuplex Ultra 360 30° 20G x 100 mm	4892510NR-01	FTR2381	25
NRFit Stimuplex Ultra 360 30° 20G x 150 mm	4892515NR-01	FTR2382	
NRFit Stimuplex Ultra 360 30° 22G x 35 mm	4892503NR-01	FTR2383	
NRFit Stimuplex Ultra 360 30° 22G x 50 mm	4892505NR-01	FTR2384	
NRFit Stimuplex Ultra 360 30° 22G x 80 mm	4892508NR-01	FTR2385	
NRFit Ultraplex 360 30° 20G x 100 mm	4892610NR-01	FTR2456	
NRFit Ultraplex 360 30° 20G x 150 mm	4892615NR-01	FTR2457	
NRFit Ultraplex 360 30° 22G x 35 mm	4892603NR-01	FTR2458	
NRFit Ultraplex 360 30° 22G x 50 mm	4892605NR-01	FTR2459	
NRFit Ultraplex 360 30° 22G x 80 mm	4892608NR-01	FTR2460	

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT

5. PERIPHERAL NERVE BLOCKS

SINGLE SHOT

HNS 12 NERVE STIMULATOR



STIMUPLEX® HNS 12 WITH SENSE

The enlarged full graphic LC-display gives you all necessary information at a glance. High contrast and a wide viewing angle ensure excellent readability.

- Stimulus amplitude in mA (large digits)
- Stimulus duration in ms
- Stimulus frequency in Hz
- Load impedance in kΩ
- Charge delivered (in nC, if activated)

BENEFITS


- Now with SENSE technology for nerve stimulation
- Advanced ergonomic shape and large display
- Digital ratchet dial for precise and tactile current adjustment
- Quick access buttons allow direct switching of key parameter values during procedure
- Alarm screens and acoustic signals inform the user of any relevant deviations

Product Description	B. Braun Code	NHS SC Code	UOM
HNS 12 Carry Case	4892095	Not Listed	
HNS 12 Electrode Cable	4892070	Not Listed	
HNS 12 Nerve Stimulator	4892098	Not Listed	1
HNS 12 Stimuplex Pen Nerve Mapping Device	4892099	Not Listed	
HNS 12 Test Resistor	2892096	Not Listed	

(Available directly through B. Braun Medical Ltd)

5. PERIPHERAL NERVE BLOCKS


CONTINUOUS NERVE BLOCKS

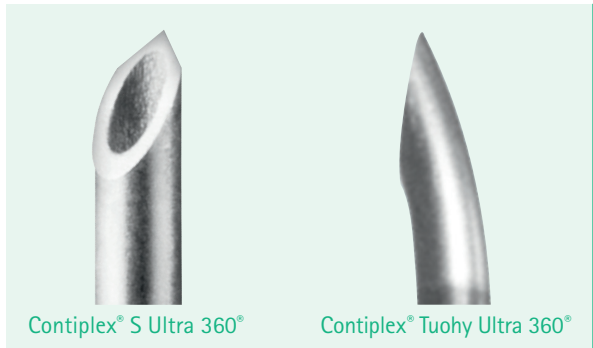

CONTIPLEX® CATHETER THROUGH NEEDLE

- NRFit™ Contiplex® S (20 degree) or Tuohy Needle
- NRFit™ Perifix® One Catheter and Thread Assist Guide
- NRFit™ Catheter Connector
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- NRFit™ Omnifix® 5 ml Lock Syringe

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Contiplex S Ultra 360 18G x 50 mm	4898650NR-27	FTR2762	
NRFit Contiplex S Ultra 360 18G x 100 mm	4898610NR-27	FTR2763	
NRFit Contiplex S Ultra 360 18G x 150 mm	4898615NR-27	FTR2764	
NRFit Contiplex Tuohy Ultra 360 18G x 40 mm	4898704NR-01	FTR2758	10
NRFit Contiplex Tuohy Ultra 360 18G x 50 mm	4898705NR-01	FTR2759	
NRFit Contiplex Tuohy Ultra 360 18G x 100 mm	4898710NR-01	FTR2760	
NRFit Contiplex Tuohy Ultra 360 18G x 150 mm	4898715NR-01	FTR2761	

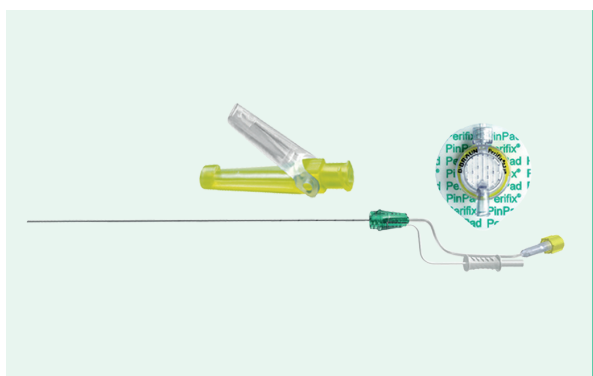

CATHETER OVER NEEDLE

- NRFit™ Contiplex® C Needle Catheter
- 30 degree or 15 degree Bevel
- NRFit™ Catheter Connector
- NRFit™ Extension
- NRFit™ Perifix® Filter 0.2µm
- Perifix® Pin Pad
- Nerve Block Sticker

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Contiplex C 25G x 190 mm 15°	4898115NR	FTR2765	
NRFit Contiplex C 25G x 190 mm 30°	4898130NR	FTR2782	5

WHAT IS NRFIT?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFIT™ ACCESSORIES

TRAINING AND SUPPORT



PNB PACK

NRFit™ CONTENTS

- 1 x NRFit™ 18G Filter Needle
- 1 x NRFit™ Ultraplex® PNB 22G x 80 cm
- 2 x NRFit™ Omnifix® 20 ml Slip

LUER CONTENTS

- 1 x Luer 18G Filter Needle
- 1 x Luer Safety Hypodermic 25G x 1.59 cm
- 1 x Luer Omnifix® 5 ml Slip

ACCESSORIES

- 5 x Swabs 7.5 cm x 7.5 cm
- 1 x Yellow Waste Bag with Adhesive Strip
- 1 x ChloraPrep 2% Chlorhexidine, 1 ml
- 1 x Blue Wrapping

Product Description

B. Braun Code

NHS SC Code

UOM

Fascia Iliaca Block Pack

BRA44100

EVS85633

40



ULTRASOUND BLOCK PACK

NRFit™ CONTENTS

- 1 x NRFit™ Needle, Blunt Fill, 18G, 50 mm (2") Yellow, Filter, Drawing Up
- 1 x NRFit™ Needle, PNB, 22G, 80 mm (3.15") 30 degree Bevel, Ultraplex® 360
- 2 x NRFit™ Syringe 20 ml Slip

LUER CONTENTS

- 1 x Needle, Blunt Fill, 18G with 5micron Filter
- 1 x Needle, Hypodermic, 25G, 0.625", Orange, Safety
- 1 x Syringe 5 ml, Luer Slip

ACCESSORIES

- 1 x Chloroprep 2% Chlorhexadine, 10.5ml Tinted
- 1 x Dressing, Adhesive
- 1 x 5 Swab, 7.5 cm x 7.5 cm
- 1 x Yellow Waste Bag with Adhesive Strip, Folded to 7.5 cm x 12 cm

Product Description

B. Braun Code

NHS SC Code


UOM

Ultrasound Block Pack

BRA45700

TBC

25



NRFIT ACCESSORIES

7. NRFit™ ACCESSORIES



NRFit™ OMNIFIX SYRINGES

- High contrast black graduation
- Plunger backstop helps to prevent unintentional plunger withdrawal
- 1 ml, 3 ml, 5 ml, 10 ml, 20 ml, 50 ml
- Slip or Lock



NRFit™ SYRINGE CAPS

- Used for storage and/or transportation of NRFit™ Syringes

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Syringe Cap	4517501N	FWC561	25

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Omnifix-F 1 ml Slip Syringe	9164001	FWC495	100
NRFit Omnifix 3 ml Lock Syringe	4617003	FWC497	
NRFit Omnifix 3 ml Slip Syringe	4616003	FWC494	
NRFit Omnifix 5 ml Lock Syringe	4617006	FWC503	
NRFit Omnifix 5 ml Slip Syringe	4616005	FWC499	
NRFit Omnifix 10 ml Lock Syringe	4617014	FWC509	
NRFit Omnifix 10 ml Slip Syringe	4616014	FWC505	
NRFit Omnifix 20 ml Lock Syringe	4617021	FWC513	
NRFit Omnifix 20 ml Slip Syringe	4616021	FTR2581	
NRFit Omnifix 50 ml Lock Syringe	4617508F-01	FWC535	



NRFit™ STERICAN 5µM FILTER NEEDLE

- Filters particles during the injection or the withdrawal of medication

NRFit™ STERICAN 5µM FILTER STRAW

- 4.5 cm or 10 cm in length

Product Description	B. Braun Code	NHS SC Code	UOM
NRFit Filter Needle 5µm	4550404N	FTR2567	50
NRFit Filter Straw 10 cm 5µm	4550200N	FTR2597	
NRFit Filter Straw 4.5 cm 5µm	4550250N	FTR2596	



WHAT IS NRFit™?

SPINAL NEEDLES AND PACKS

COMBINED SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE BLOCKS

NRFit™ ACCESSORIES


TRAINING AND SUPPORT



SCAN QR CODE FOR MORE INFORMATION ON NRFit™ SYRINGES




8. TRAINING AND SUPPORT

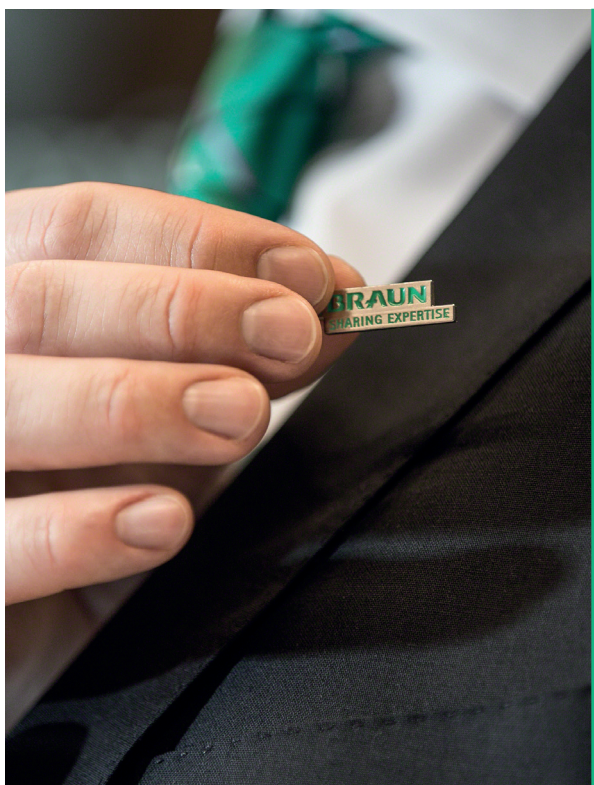

We offer comprehensive and educational support during implementation including:

- Multiple B. Braun staff covering all clinical areas
- Removal of luer stock
- Stocking of NRFit™ products into trolleys
- Hands-on training sessions
- Awareness posters (see example)
- Training posters (see example)

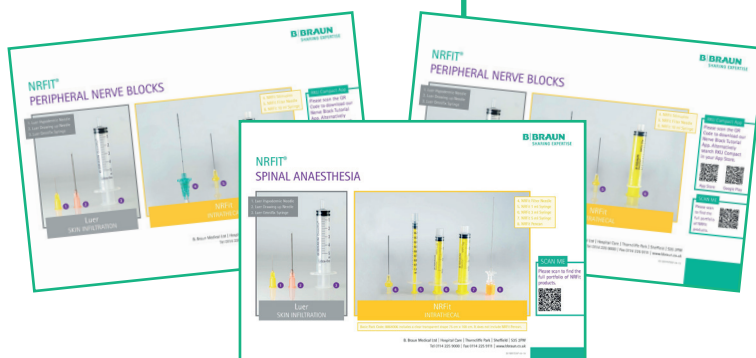

OUR MISSION

We develop
EFFECTIVE
solutions and
guiding standards for the
HEALTHCARE
system through
CONSTRUCTIVE
dialogue with our
customers and partners.

WHAT IS
NRFT™?

SPINAL NEEDLES AND
PACKS

COMBINED
SPINAL EPIDURAL

EPIDURALS AND PACKS

PERIPHERAL NERVE
BLOCKS

NRFT™
ACCESSORIES

TRAINING AND
SUPPORT

References:

1. <https://www.england.nhs.uk/wp-content/uploads/2014/02/non-Luer-spinal-suppl.pdf>
2. Kåsine, T., Romundstad, L., Rosseland, L., Ullensvang, K., Fagerland, M., Hol, P., Kessler, P. and Sauter, A., 2019. Needle tip tracking for ultrasound guided peripheral nerve block procedures—An observer blinded, randomised, controlled, crossover study on a phantom model. *Acta Anaesthesiologica Scandinavica*, 63(8), pp.1055-1062.
3. Kåsine, T., Romundstad, L., Rosseland, L., Ullensvang, K., Fagerland, M., Kessler, P., Bjørnå, E. and Sauter, A., 2019. The effect of needle tip tracking on procedural time of ultrasound guided lumbar plexus block: a randomised controlled trial. *Anaesthesia*, 75(1), pp.72-79.